


# Camera/Footswitch Covers

## Camera Covers



Mediseal® camera covers are designed to protect a wide variety of camera sizes.

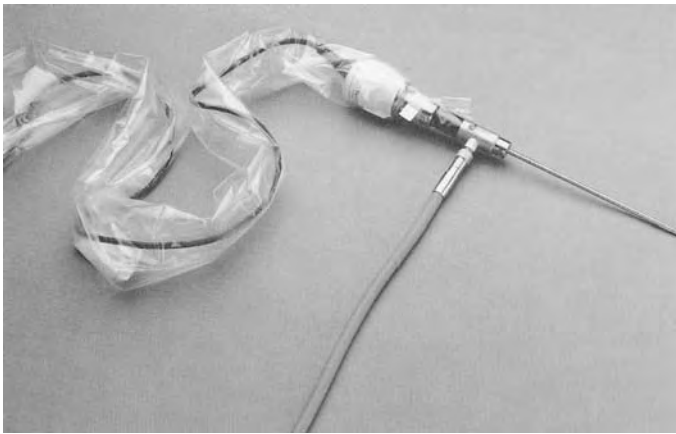

**#IG-CK-406** Camera Cover  
7" Proximal, 4" Camera, 72" Long,  
2 3/4" Tip, Supplied Sterile  
10/Bx

While using the Mediseal® camera covers, you will find that the embossed plastic will give you a secure grip while allowing you to have better and easier control of your camera and laser.

Custom-made for cameras, Mediseal® covers are wider at the proximal end and are telescopically folded for easy camera insertion. The stepped-down distal tip, with a shoulder smaller than the camera body, prevents the camera from passing through the cover and into the sterile field.

With convenient adhesive strips and twist ties, the Mediseal® cover can be secured to the scope easily maintaining the sterile field. The Mediseal® camera covers will allow you to work with confidence while knowing that contaminants are not entering the O.R. sterile environment.

## Footswitch Covers

The Mediseal® drawstring makes the footswitch covers easy and fast to secure, saving time and reducing worry. Just slip the footswitch into the cover, tighten the drawstring, loop it around the neck of the cover and tie with a simple knot.

The clear, durable Mediseal® plastic cover makes footswitches easy to see and even easier to use. The covers are proven to withstand the stress of repeated foot action.



**#IG-FC-904** Footswitch Cover  
14" Wide, 24" Long, 4" Gusset  
Non-Sterile  
20/Bx

Your equipment works day in and day out, reducing footswitch failure and the resulting costly delays and frustration in O.R. scheduling.

